


ELECTROSURGICAL ACCESSORIES

Microsurgical Needle

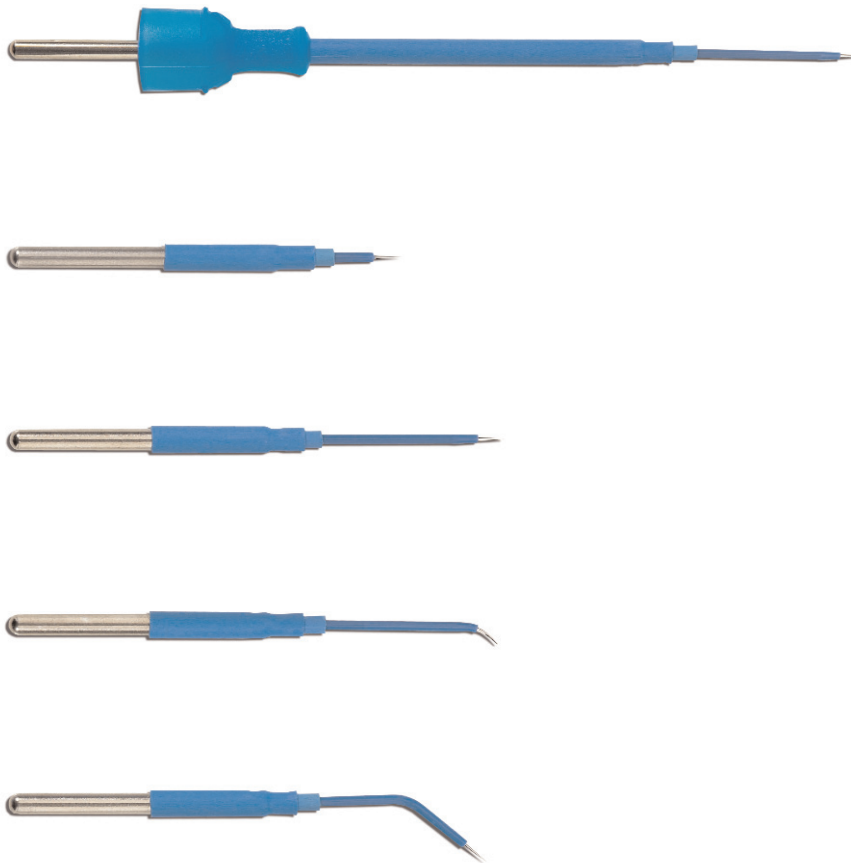

Introducing microsurgical tungsten needles from Valleylab. These electrodes offer improved cutting and coagulation performance, while minimizing thermal spread. These microsurgical needles offer surgeons the precision they need for fine dissection in a wide variety of procedures.

PRECISION

The needles' heat-resistant tungsten tips maintain sharpness throughout the procedure, enabling fine, precise tissue dissection at very low power settings.

ENHANCED PERFORMANCE

These microsurgical needles offer greater precision at lower power settings than traditional stainless steel needles in both Cut and Coag modes, virtually eliminating electrode drag and reducing thermal damage.

IMPROVED SAFETY

The needles' protective high temperature insulation extends to within 3 mm of the electrode tip, guarding against accidental tissue contact and burns. The unique Safety Sleeve™ feature of the E1465B reduces the risk of alternate site burns, especially in procedures like tonsillectomies.

VERSATILITY

These needles are compatible with most electrosurgical pencils and come in a selection of five lengths and angles, allowing the surgeon to choose the electrode appropriate for each procedure.




Microsurgical Needle


Order Information

CATALOG NUMBER	DESCRIPTION	ORDER QUANTITY
----------------	-------------	----------------


Sterile, single use electrodes
shown actual size

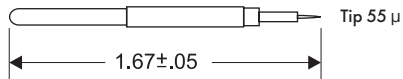


E1650 2 cm straight microsurgical tungsten needle 10/case

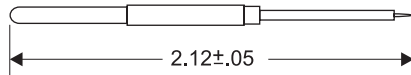


E1651 3 cm straight microsurgical tungsten needle 10/case

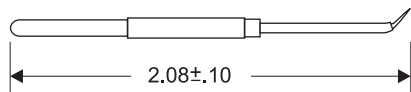

E1652 3 cm microsurgical tungsten needle, 45 degree angle, 3 mm from tip 10/case

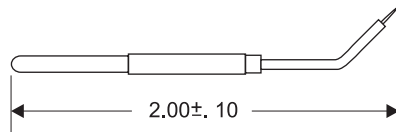

E1653 3 cm microsurgical tungsten needle, 45 degree angle, 1 cm from tip 10/case

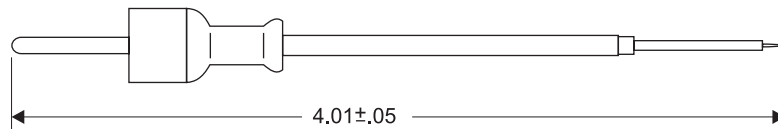

E1654B 9 cm (3.5 in.) straight microsurgical tungsten needle with Safety Sleeve™
for use only with Valleyslab handswitching electro-surgical pencils 10/case

Maximum power of 25 watts for use with microsurgical needle

Valleyslab
a division of Tyco Healthcare Group LP
5920 Longbow Drive
Boulder, CO 80301-3299 USA
800 255 8522

