


ELECTROSURGICAL ACCESSORIES

Single Use Arthroscopic Electrodes

Performance, Safety, and Quality

LIGASURE™ VESSEL SEALING SYSTEM
ELECTROSURGICAL GENERATORS
ELECTROSURGICAL ACCESSORIES
LAPAROSCOPIC INSTRUMENTS
SMOKE EVACUATION PRODUCTS
ARGON-ENHANCED ELECTROSURGERY
COOL-TIP™ RF ABLATION SYSTEM
RF PAIN MANAGEMENT SYSTEM
EDUCATION & TRAINING

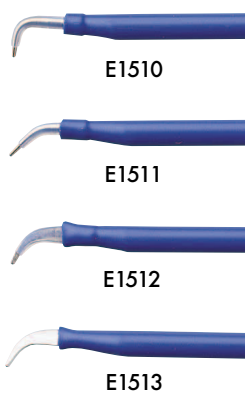

Arthroscopic electrodes from Valleylab offer enhanced performance and safety features. Available in four tip configurations and sold separately or packaged with a Valleylab electro-surgical pencil, arthroscopic electrodes bring Valleylab performance, quality and reliability to arthroscopic procedures.

PRECISE ENERGY DELIVERY

Transparent insulation extends close to the electrode tip, improving the surgeon's visualization, allowing the use of lower power settings and performing well in the conductive irrigating solutions used frequently in arthroscopy.

ENHANCED SAFETY

The electrode cap prevents fluids from entering the electro-surgical pencil, stabilizes the electrode, and ensures a safe connection between the electrode and pencil.

COMPATIBLE WITH THE VALLEYLAB ELECTROSURGERY SYSTEM

Use arthroscopic electrodes with the Valleylab Force FX™ generator with Instant Response™ technology to achieve a smooth precise cut, consistent clinical effect through all tissue types, and decreased peripheral tissue damage at lower power settings. Valleylab REM™ patient return electrodes virtually eliminate the risk of patient burns at the return electrode site, and Valleylab electro-surgical pencils are renowned for their quality and reliability.



Single Use Arthroscopic Electrodes

Technical Specifications

Total Length:	21.6 cm (8.5 in)	
Active Length:		
E1510	0.19 cm (0.075 in.)	
E1511	0.19 cm (0.075 in.)	
E1512	0.19 cm (0.075 in.)	
E1513	0.7 in (0.275 in.)	
Diameter of Shaft:	2.9 mm 0.115 in	
Diameter (including tip):		
	Min	Max
E1510	4.45 mm (0.175 in.)	4.95 cm (0.195 in.)
E1511	3.81 mm (0.150 in.)	4.57 mm (0.180 in.)
E1512	4.06 mm (0.160 in.)	4.83 mm (0.190 in.)
E1513	4.06 mm (0.160 in.)	4.70 mm (0.185 in.)
Insulation:		
Shaft:	Polyolefin, Blue	
Tip:	PTFE, Natural	

Compatible with: All Valleylab handswitching electro-surgical pencils.

Order Information

CATALOG NUMBER	DESCRIPTION	ORDER QUANTITY
Rockerswitch Pencil with Electrode		
E2510H	Arthroscopic Hook Electrode, insulated	10/case
E2511H	Arthroscopic Hook Electrode, 45° Angle, insulated	10/case
E2512H	Arthroscopic Angled Blade Electrode, insulated	10/case
Electrodes		
E1510	Arthroscopic Hook Electrode, insulated	10/case
E1511	Arthroscopic Hook Electrode, 45° Angle, insulated	10/case
E1512	Arthroscopic Angled Blade Electrode, insulated	10/case
E1513	Arthroscopic Angled Blade Electrode, non-insulated	10/case



Valleylab
a division of Tyco Healthcare Group LP
5920 Longbow Drive
Boulder, CO 80301-3299 USA
Valleylab Direct: 800 255 8522
To Order Product: 800 722 8772

tyco / Healthcare